

Department of Physics and Astronomy
University of Heidelberg

Bachelor Thesis in Physics
submitted by

Robin Eberhard
born in Aalen, Germany

handed in on
March 29, 2016

Characterization of a multispecies imaging system

This Bachelor Thesis has been carried out by Robin Eberhard at the
Physikalisches Institut in Heidelberg
under the supervision of
Prof. Dr. Matthias Weidemüller

Characterization of a multispecies imaging system

Robin Eberhard

Abstract This work describes the characterization and implementation of a new camera system for imaging of atomic density distributions in a multispecies ultracold atom experiment. The high quantum efficiency of the camera and its new readout mode, which allows for faster acquisition, are perfectly suited for scientific imaging. In this work, the noise sources that are induced from the camera itself, being the readout and dark noise, are described. To understand these, the setup of a typical CCD-camera is explained. At last, the imaging setup, including the magnification of an atomic cloud was tested on an ideal Fermi gas. The high resolution of the system and the noise reduction shows a precise distribution, which is explained and the degeneracy parameter is then extracted to be around $T/T_F =$

temperature...

Zusammenfassung Lorem ipsum dolor sit amet, consectetur adipiscing elit. Ut purus elit, vestibulum ut, placerat ac, adipiscing vitae, felis. Curabitur dictum gravida mauris. Nam arcu libero, nonummy eget, consectetur id, vulputate a, magna. Donec vehicula augue eu neque. Pellentesque habitant morbi tristique senectus et netus et malesuada fames ac turpis egestas. Mauris ut leo. Cras viverra metus rhoncus sem. Nulla et lectus vestibulum urna fringilla ultrices. Phasellus eu tellus sit amet tortor gravida placerat. Integer sapien est, iaculis in, pretium quis, viverra ac, nunc. Praesent eget sem vel leo ultrices bibendum. Aenean faucibus. Morbi dolor nulla, malesuada eu, pulvinar at, mollis ac, nulla. Curabitur auctor semper nulla. Donec varius orci eget risus. Duis nibh mi, congue eu, accumsan eleifend, sagittis quis, diam. Duis eget orci sit amet orci dignissim rutrum.

Contents

1. Introduction	1
2. Setup for high resolution imaging	3
2.1. Experimental requirements	3
2.2. Camera for double species imaging	4
2.2.1. Basics of CCD cameras	4
2.2.2. Comparison with the previous setup	6
2.2.3. Dark current	7
2.2.4. Readout noise	8
2.2.5. Quantum efficiency	10
2.3. Mechanical shutter	11
2.3.1. Electronic and mechanical setup	11
2.3.2. Dynamical properties	13
2.4. Mask for the CCD sensor	14
2.4.1. Fast kinetics mode	15
2.4.2. Frequency response of a slit	18
2.4.3. Optimization of the masking setup	21
3. Thermometry of an ultracold ideal fermi gas	23
3.1. Absorption imaging	23
3.2. Density distributions of ideal Fermi gases	24
3.3. Finding properties of the Fermi gas	26
4. Conclusion and outlook	29
A. Acquisition sequence	31
B. Design of the shutter	33
C. Shutter circuit	35
D. Setup of the custom slit	37

1. Introduction

Conducting experiments, measuring attributes and taking data are quantities that go hand in hand in scientific fields. Already in the early 1900s, imaging systems were chosen as data acquisition in astronomy. A complete optical system allowed for magnification of planets, our sun or nebulae. Already then, they were able to store their acquisitions on photo plates [1], while nowadays so-called charge-coupled devices (CCD) or complementary metal-oxide-semiconductors (CMOS) are used as they can take the data digitally.

This also found application in the field of ultracold quantum physics — an experimental field that is around since the 1970s. Hereby, atoms are prepared in traps in order to reduce their velocities. The atoms can then be detected by illuminating them with an imaging beam, a laser beam that is resonant with the atoms. The atoms will cast a shadow, which can be focused onto the chip, where the result is digitalized and can be analyzed. Interesting attributes of atomic clouds are usually density distributions and atom numbers. For the latter, a high resolution imaging system is not needed, since the optical density is extracted from the whole chip and not individual pixels. Therefore if the magnification of the system is known, one can directly conclude the atom number.

In order to directly see interactions between atoms and their distributions they yield, a high resolution system is commonly used. Therefore a system of lenses needs to be carefully set up in order to gain maximum resolution. A high number of pixels on the chip then allows to see detailed information of the cloud. Especially for small number of atoms it is additionally important to have small noise in the camera as well as a high sensitivity to photons, being the quantum efficiency.

Interesting features were already observed for example in fermionic superfluids, where a gas of spin imbalanced Fermions were prepared. The cloud was ”stirred” using a laser beam, revealing small vortices inside the cloud which were visible as holes.

Our experiment uses a mixture of ^6Li and ^{133}Cs . The new imaging system is particularly well suited for the study of the Efimov effect [2–4] and different polarons (Fermi-polaron, Bose-

polaron, Bubble-polaron). In case for the Bubble polaron, the resolution of the imaging system will allow to find a hole in the density distribution of the atoms.

This thesis builds upon the optical system, that made the high resolution possible and will introduce a new camera into the experiment, which has new features, that allow for noise reduction and an improvement in acquisition timings. At first, the imaging setup is explained in Chapter 2.1. Chapter 2.2 then introduces the camera used and will highlight its qualities and improvements to the system2.2.2. In order to fully understand noise accumulation on the chip, the setup of a CCD is then explained in Chapter 2.2.1, followed by different noise sources being the dark current2.2.3, which is a temperature dependent noise that will accumulate over time and the readout noise2.2.4, which adds noise as the image is being read out. Subsection 2.2.5 will then give an overview about the quantum efficiency and how it affects the imaging quality.

In Chapter 3 a measurement was concluded in order to test the camera. An ideal Fermi gas was prepared, which shows a unique distribution different from Bose-Einstein condensates, from which the temperature and Fermi temperature were concluded.

2. Setup for high resolution imaging

In order to understand the nature of atoms in more detail, experiments are conducted where they are trapped and cooled. Under these circumstances, interactions between atoms can be tuned using Feshbach resonances. This allows to build up for example Bose-Einstein condensates (BEC) or Cooper pairs. The complete setup has been explained in more detail before [5]. The following sections will focus on the imaging in order to extract these attributes and the measures it takes to reduce noise sources on the camera.

2.1. Experimental requirements

In this double species experiment, fermionic ${}^6\text{Li}$ and bosonic ${}^{133}\text{Cs}$ are trapped and cooled inside a vacuum chamber. The atoms are emitted from an oven into a Zeeman slower, which transfers them into a magneto-optical trap. The MOT will further cool the atoms up to the doppler limit (around $140\text{ }\mu\text{K}$). To reach lower temperatures, the atoms will be trapped inside an optical dipole trap. Forced evaporative cooling will release the fastest atoms, until temperatures of a few nano Kelvin are reached.

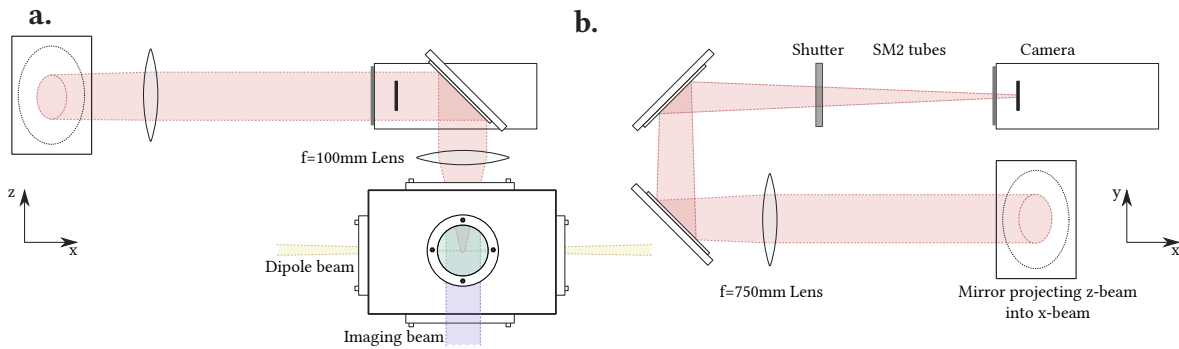


Figure 2.1.: **Imaging path.** The camera is mounted above the vacuum chamber as visible in **a.** (gravity would be in negative z direction). The image of the atoms is collimated with the first lens and passing the optics in **b.**, until it is refocused in the chip on the camera.

In order to image the atoms, an imaging beam is pointing along the z -axis according to Figure 2.1. An achromatic doublet lens focuses the imaging beam onto the CCD camera, which is located on top of the vacuum chamber.

In order to measure even at low atom numbers ($n = 1000$), the setup was refined allowing for high resolutions. The first lens, which collimates the image of the atoms has a low focal length ($f=100$ mm). The image is refocused into the camera by the second lens. Since the imaging uses two separate frequencies to image both Lithium and Caesium, the chromatic shift introduced by the lenses will result in two different focal points for the images. This can be compensated by a high focal length, which therefore allows to have the camera in the Rayleigh range of both laser beams, when they are assumed to be gaussian.

Having the first lens as close as possible to the atoms then allows for the highest resolution. With the high resolution of the camera it is then possible to take detailed images of atomic clouds.

2.2. Camera for double species imaging

The camera used in order to find atoms, the Andor iKon M [6], is a charge-coupled device (CCD). Using this camera, noise sources such as Readout noise and Dark noise can be reduced, by cooling the chip, which are explained in the following sections. The camera also ships with a new acquisition mode, where it is now possible to take two consecutive images of the same atom cloud in a time of flight (ToF) measurement.

2.2.1. Basics of CCD cameras

A camera operates by means of converting photons first into electrons then into voltage, which is finally read out as data [7] [8]. Each conversion process can add noise to the final image, which needs to be minimized in order to acquire accurate data.

The photons are collected on an array of semiconductor photo diodes, called the pixels, where ideally the spacing between the pixels is zero to gain maximum accuracy. The resolution is then dependent on the pixel size, which is usually between $10\text{ }\mu\text{m}$ and $20\text{ }\mu\text{m}$ per pixel for scientific cameras. Bigger pixels mean higher photon sensitivity but usually lower resolution.

To create a digital image, the charges from the pixels have to be shifted one-by-one into the analog to digital converter (ADC). This is done by vertically shifting them into the

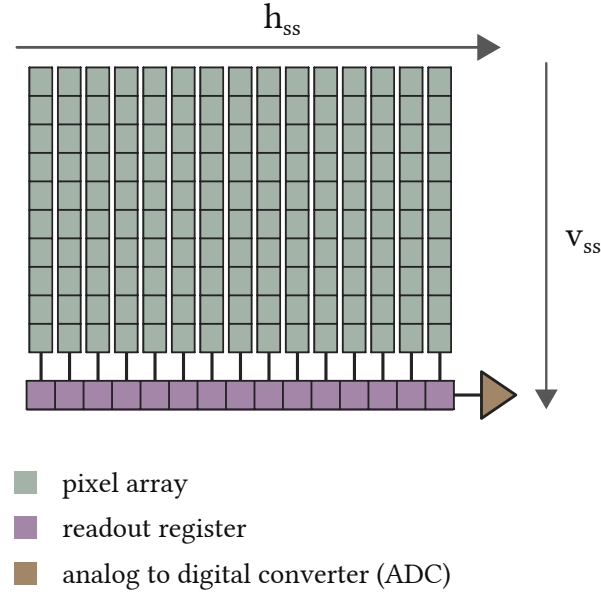


Figure 2.2.: **Schematic design of a CCD array and its readout.** The pixels are arranged in the pixel array. During readout they are shifted into the readout register and then to the side into the analog to digital converter (ADC).

readout register and then horizontally into the ADC, where the charges are multiplied and converted to digital data. This is illustrated in Figure 2.2.

The shifting is done by storing the charges after collecting them, where each storage can be seen as an electronic potential well. In order to shift the charges and prevent overlapping, three spatially separated regions with potentials U_1 , U_2 and U_3 are required. Figure 2.3 indicates the systematics behind the shifting.

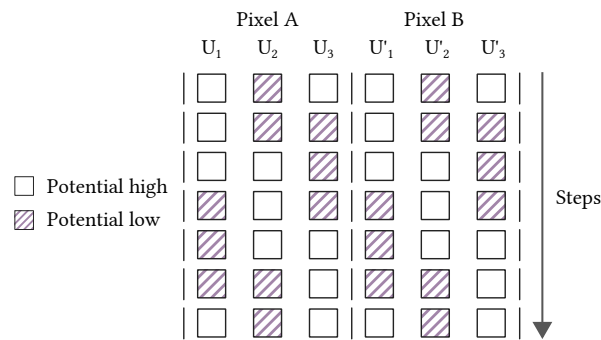


Figure 2.3.: **Shifting charges in a CCD detector.** To shift charges from Pixel A to Pixel B, the three potentials at the regions U_1 , U_2 and U_3 in each pixel have to be set accordingly to allow the charge flow without overlapping each other. Each row represents a single step.

Moving electrons to the next region is a three-step process. The charges are first only

present in the region U_2 , while the potentials U_1 and U_3 are kept high. They are then distributed across U_2 and U_3 by setting them both low. At last U_2 will be set high such that the charge is now fully in U_3 . This is repeated until the charge has been finally shifted from pixel A to B, which in total takes seven steps.

2.2.2. Comparison with the previous setup

For the imaging of small atomic clouds, it is very important to have cameras with minimal noise and maximal readout speed. The new setup improves both. The dark noise, which is a temperature dependant effect, where counts accumulate over time, can be significantly reduced by cooling the chip down to less than -70°C . A new readout mode, called fast kinetics, makes it possible to acquire all images before reading out. This significantly improves the speed at which images can be taken.

The readout speed is highly important in our setup. Since absorption imaging is the technique of choice to measure atom attributes, three images need to be taken in each sequence, being the absorption, division and background image. The old setup used a Guppy-38B [9] camera, which has a frame rate of 30fps. This meant the acquisition was finished after 100 ms. The new Andor camera, on the other hand, can take images quickly without the need to read out in-between. At the fastest shift speed, the acquisition is finished after 1.632 ms, improving the speed by a factor of more than 60.

The quantum efficiency is explained in Chapter 2.2.5 and is also highly important, as it describes how many photons are detected on a camera and is directly connected to the sensitivity. A higher quantum efficiency therefore means better results. They compare as $QE_{Li} = 79\%$, $QE_{Cs} = 77.3\%$ for the Andor and $QE_{Li} = 35\%$, $QE_{Cs} = 10\%$ for the Guppy camera, therefore the sensitivity of the new camera is significantly higher.

When comparing the resolution, the chip size also has to be considered. Since higher resolutions seem to be preferable at first, it also means that for the same pixel sizes, the photon sensitivities will decrease. The pixel sizes are $8.4\ \mu\text{m} * 9.8\ \mu\text{m}$ in the old setup and $13\ \mu\text{m} * 13\ \mu\text{m}$ in the new setup, while the resolutions are 768×492 and 1024×1024 respectively. This means that the new setup allows larger magnification.

The Guppy camera is a lot smaller than the Andor ($48.2\ \text{mm} \times 30\ \text{mm} \times 30\ \text{mm}$ vs $204.2\ \text{mm} \times 105\ \text{mm} \times 107\ \text{mm}$), therefore making an implementation on a full experimental table easier. To implement the complex imaging system surrounding the Andor camera, a lot of preparation was made in the thesis of Carmen Renner [10].

In spite of its large size, the new camera offers the ability to image both Lithium and Caesium species at once, while two Guppy cameras were needed beforehand, which also meant placing them on different imaging axes.

2.2.3. Dark current

A common noise source that is apparent in all CCD cameras, is the so-called dark current. It originates from the thermal excitation of electrons in individual pixels. Since they are made of semiconductors, once in a while, an electron can pass the potential between valence and conduction band simply due to their thermal energy. Thus, excess electrons accumulate, which contribute to the background signal and introduce additional noise [11].

The dark current has a strong temperature dependence [12]

$$I_{dark}(T) \propto T^{\frac{3}{2}} \exp(-E_g/2k_B T), \quad (2.1)$$

where E_g is the band gap, that separates the valence from the conduction band in the semiconductor, T the temperature and k_B the Boltzmann constant. Therefore in order to reduce dark current noise, the temperature of the chip can be reduced, which decreases the thermal energy of the electrons.

This reduction of dark current with lower temperature has been measured and verified in Figure 2.4, using the built-in peltier element in the Andor camera to control the temperature of the CCD chip. For a long exposure time, dark current accumulates on all pixels. The counts are measured for several temperature settings and the counts C are converted to electrons per pixel per second (I_i) [12] for each pixel i as

$$I_i = \frac{C}{G t_{exp}}, \quad (2.2)$$

with the Gain G of the camera, and the exposure time t_{exp} . The readout time of the pixels can be neglected, as it is significantly lower than the exposure time (see Chapter 2.4.1). The data should then follow the theory in Equation 2.1, although deviations are visible in the low temperature regime, which result from the fan not being able to divert the heat from the chip.

As it is important to minimize the noise source, the chip is cooled to the lowest possible temperature. This reduces the possibility of thermal electrons creating counts in the detector, so that even at small atom numbers, they are still distinguishable from the noise.

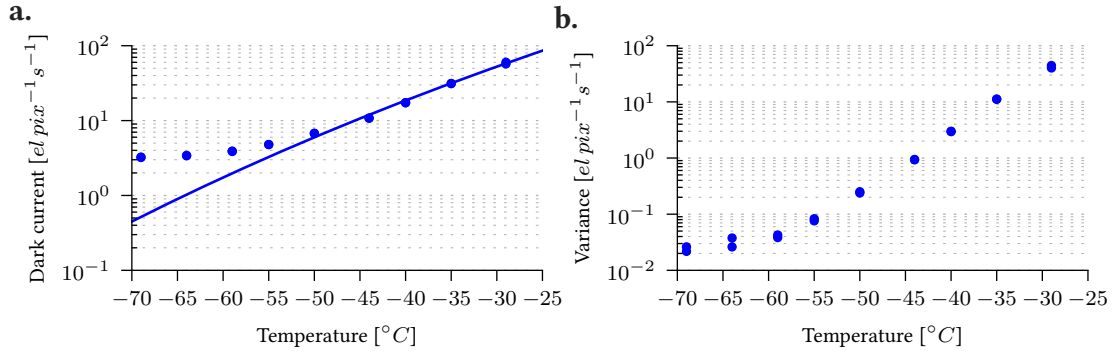


Figure 2.4.: **Dark noise accumulation on the chip.** For a long exposure time of 100 s, the dark current was measured for several temperatures. Gain in this measurement was minimal (0.215 [10]). Equation 2.1 was fitted, leaving the amplitude A and band gap E_g free, which resulted in $A = 1.1 * 10^8 \text{ W m}^{-1} \text{ K}^{-1}$ and $E_g = 0.96 \text{ eV}$.

2.2.4. Readout noise

As described in Chapter 2.2.1, pixels are shifted in order to be read out by the ADC. Moving charges from pixel to pixel causes noise that accumulates over each iteration. Together with the noise which the ADC introduces, this is then called the readout noise. The shifting noise together with the dark noise is visible as a gradient (Figure 2.5) since each shift adds new charges due to excitations in the semiconductors. In order to characterize this, one can take the variance of the image, which should be zero for no noise.

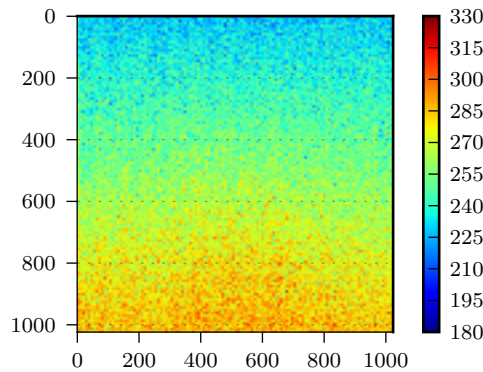


Figure 2.5.: **Gradient on the chip during readout.** Dark current accumulates on the chip over time. During readout, this is visible as a gradient, since the pixels that have not been read out are still affected by this noise. As the chip temperature is being lowered, the gradient will also degrade.

In Figure 2.6, the readout speeds were probed for their noise. The measurement was carried

out by taking images at low exposure time and low temperature, so that the main noise source would be the readout and not the dark current. The variance of an image then indicates the noise accumulation on the chip and it has been shown, that slower readout speeds accumulate less noise, although it is only dependent on the horizontal, not on the vertical shift speed.

The actual readout of the chip can be done at slow speeds, so that the noise from horizontal shifting is low. As the vertical shifts do not add significant noise, it can safely be set to fast speeds, which is also favorable, as ^6Li and ^{133}Cs are imaged back to back as fast as possible, so that the cloud dynamics are still approximately the same.

Shifting charges means increasing and decreasing potential wells, so that they can move from one into the other. If this is done fast enough (non-adiabatic), electrons in the pixel are more likely to get excited from the valence to the conduction band from the semiconductor, therefore adding noise. As a consequence of these results, slower readout speeds are used in order to not add significantly more noise.

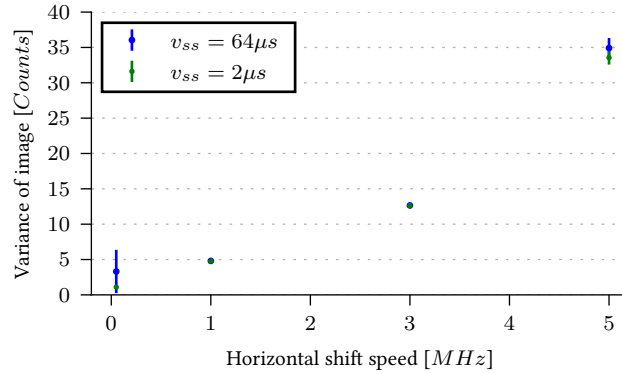


Figure 2.6.: **Readout noise.** The pixels are shifted row-wise into the readout register, depending on the vertical shift speed (v_{ss}) and then moved pixel-by-pixel with the horizontal shift speed into the analog to digital converter. Since noise reduction is important, minimal horizontal shift speeds will be used, while the vertical shift speed does not seem to affect the variance. To make the readout the dominant noise source, temperature was set to -69°C and exposure to 10 ms. The results have been received from taking the mean and variance of multiple sets of exposures.

2.2.5. Quantum efficiency

When selecting cameras for scientific imaging, one of the attributes to look out for is quantum efficiency (QE). Not all incoming photons are converted into electrons, which adds an uncertainty and additional noise to measurements. This is characterized by quantum efficiency, which is defined as the fraction of photons that are converted into electrons. A QE of 100% would mean, that every incident photon is converted into an electron on the chip. This can be put into equations as

$$QE = \frac{N_{detected}}{N_{total}}, \quad (2.3)$$

where $N_{detected}$ describes the number of detected photons and N_{total} the number of photons implying on the detector.

The measurement of the quantum efficiency has been previously carried out by Carmen Renner [10] in her diploma thesis. In order to do so, the number of detected photons can also be expressed in terms of the energy E_{beam} of all photons reaching the camera and the energy $E_\gamma = h\nu$ of a single photon [13]:

$$N_{total} = \frac{E_{beam}}{E_\gamma} = \frac{Pt_{exp}}{E_\gamma}, \quad (2.4)$$

where P is the beam power and t_{exp} the exposure time.

Now the quantum efficiency reads

$$QE = \frac{h\nu N_{detected}}{Pt_{exp}}. \quad (2.5)$$

In order to find the QE, one would now measure the accumulated counts for several exposure times, since the frequency ν and the power P are known. For the Andor camera used in our experiment, they have been measured as 70.7 % for 671 nm and 69.9 % for 852 nm.

This camera was especially chosen for its quantum efficiency, since it is important in our case, to be able to detect most photons from Lithium and Caesium absorption.

2.3. Mechanical shutter

Due to the high quantum efficiency, the camera is very sensitive to stray light. Therefore it is necessary to cover the chip from stray light. A mechanical shutter which can be electronically controlled is therefore built in front of the light path, which prevents the camera from damage between measurements. The electronics were set up and the shutter characterized in order to find perfect timings, so that light only enters the chip, when a measurement is running.

2.3.1. Electronic and mechanical setup

The high QE of the Andor camera translates into a high sensitivity of the CCD chip to stray light. Therefore, in order to not unnecessarily illuminate the chip between measurements, a mechanical shutter was built into the optical path, which can also be seen in Figure 2.1. By opening the shutter shortly before the imaging sequence and closing it immediately afterwards, the stray light that is impinging the CCD detector can be minimized.

The shutter has five fans, as seen in Appendix B, which are mechanically guided, such that they together perform a circular motion outwards. The circular motion is best achieved for narrow fans, therefore needing more in order to close the shutter properly.

The guides are connected to each other and can be pulled outwards with a mechanical switch, that can be manually pushed or pulled. To drive the shutter electronically, a magnetic coil and a magnet are used. When the coil receives a current, the magnet, which is connected to the switch, will pull the guides, moving them outwards, therefore opening the shutter.

Although the coil can pull the magnet in, it cannot be pushed away. This is compensated using a spring, which is connected to the switch. This also means, that the current driving the coil needs to be high enough to also work against the spring.

The optimal case would be now to have fast opening and closing times, since we want to prevent illumination between measurements. This can be achieved by testing several springs and to make the pull from the magnetic coils as fast as possible, which means increasing the current to drive it.

A custom circuit as shown in Figure 2.7 is used. In order to obtain a high current to increase the pull from the magnetic coil on the magnet, a transistor is used, which is controlled by a MOSFET driver. The complete circuit also contains a voltage regulator, so that the user

does not need to know the input of the driver and can simply use a high voltage to open the shutter ($V \in [2, 8] \text{ V}$) and a low voltage ($V = 0 \text{ V}$) to close it.

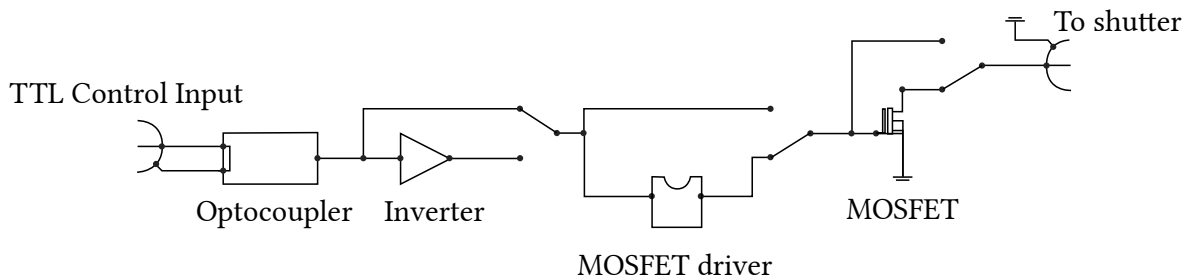


Figure 2.7.: **Electronic circuit to control the shutter.** The first element in the circuit is the optocoupler which is used to decouple parts from the circuit in order to not create loops, which would induce charges when the magnetic fields are on. The remaining parts can be skipped, by setting the jumpers. The inverter will flip the sign on the voltage. The MOSFET driver is used in order to serve the correct voltage at the transistor's gate and to help during discharge, such that the MOSFET will not overheat. The transistor will finally serve a high current of 2.5 A, that is throughput to the shutter. The jumpers in this figure are set as they are used in the experiment. The complete circuit can be found in Appendix C.

2.3.2. Dynamical properties

The shutter operates by controlling a magnetic coil, which pulls a magnet. The magnet is reverted into its original position by a spring, closing the shutter. It was discussed before, in Subsection 2.3.1, that more fans give a better approximation to a circular motion. Optimally, we would expect the shutter to open perfectly circular with linear velocity. Therefore, to optimize on this, one would look at opening and closing speeds, which should be minimal.

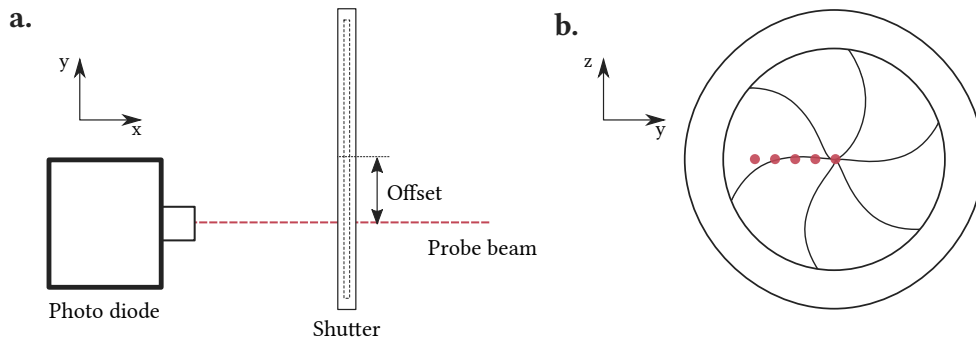


Figure 2.8.: **Probing the shutter for dynamics.** The shutter is probed at several positions using a laser. The offset in y direction was varied to find opening and closing times as a function of their offset.

To find the actual dynamics of the shutter, an experiment was set up. The shutter was probed at several positions, using a laser, to find the opening time. An example of one of many measurements is shown in Figure 2.9. The laser beam has a gaussian intensity distribution, which results in a error function on the photodiode as only partial intensity is received from the diode, when the shutter is still blocking parts of the beam.

The points have been found by pointing a laser beam at a photodiode, which was blocked by the shutter. Since the laser has a finite radius, there is a transition in the signal from the minimum to the maximum, which is due to the approximately linear opening velocity of the shutter and the gaussian intensity distribution of the laser beam.

In order to find the opening and closing times of the shutter as seen in Figure 2.10a., the offset from the centre of the shutter was varied. For larger offsets, the time from the initial trigger until a signal in the photodiode was received is longer when the shutter is opening, while it is shorter when it is closing. In Figure 2.10b. the speed is not perfectly linear. The deviations close to the centre originate from the shape of the fans. As there are five fans, there is no mirror symmetry, so that opening to the left and right is different in the end.

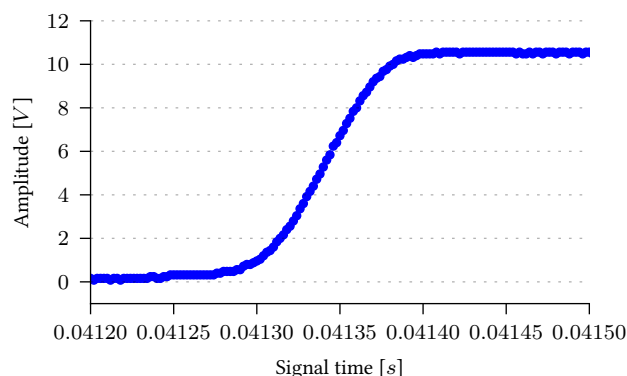


Figure 2.9.: **Shutter characterization.** The dynamics of the shutter were measured using a laser beam with a variable horizontal offset, which is fixed in this plot, and a photo diode measuring the laser intensity. In this figure, the offset is 9.3 mm from the center. An error function was fitted yielding the time until the shutter opens. The opening time is defined here as the mean value of the error function, and was found to be 41.35 ms in this case.

The results will later be used to time the triggers in the measurement, so that the shutter can stay closed as long as possible. But in order to optimize it, there were two shutters at hand, and several springs to choose from, the combination of which were all tested until the optimal timings for opening $t_{open} = 0.12$ s and for closing $t_{closing} = 0.14$ s were found.

It has been found that the opening time actually consists of three timings. The initial trigger signal runs through the electronics, until it is sent to the shutter coil. The coil then has a delay, until it starts pulling the magnet. The opening time is then found by additionally adding the time until the fans are fully open. This is also shown in Figure 2.11.

2.4. Mask for the CCD sensor

One of the key features of the Andor iKon M is the fast kinetics readout mode. This allows for fast acquisition, which is important when imaging multiple species. Using this acquisition mode, we are able to take images $500\text{ }\mu\text{s}$ apart from each other, so that it is possible to take consecutive images from the same cloud before it is too dilute. The acquisition mode is explained in these sections as well as the diffractions introduced by a slit, which is needed to mask parts of the chip.

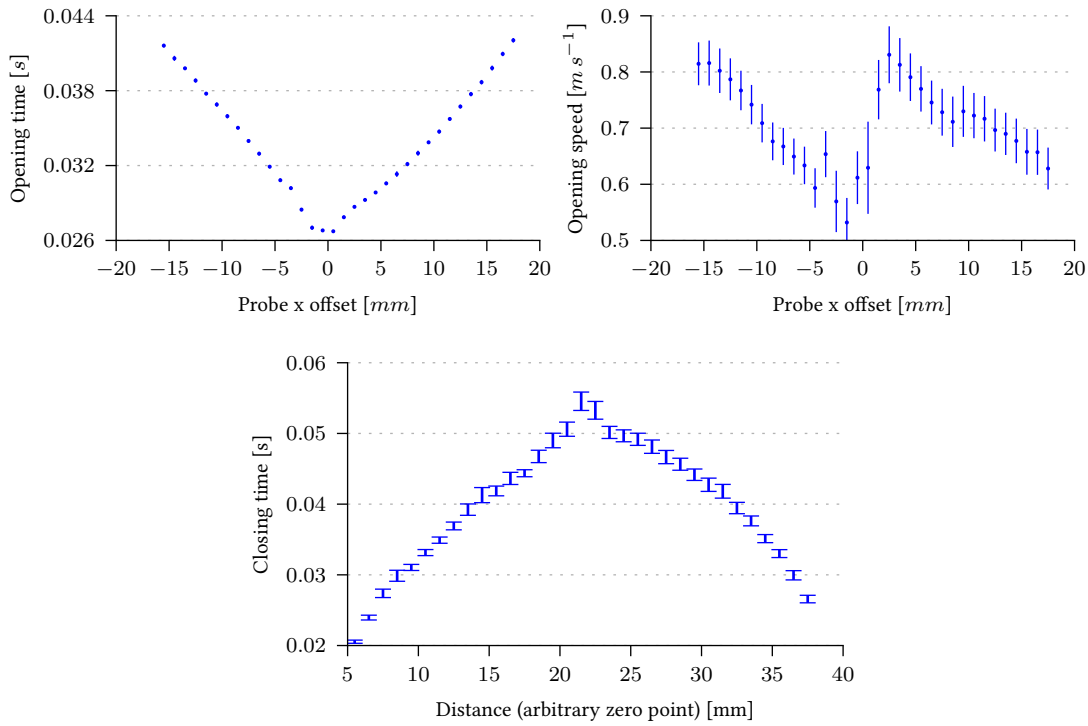


Figure 2.10.: **Sample dynamics.** The figure shows opening and closing time as well as the opening velocity for a strong spring. The spring closed the shutter quickly while it prevented it from opening fast. The velocity was measured by using the beam diameter as the distance the shutter needed to transverse. For each offset, a set of 100 images were taken and the errors found as being the variance. It is noticeable, that the opening velocity on the right side is faster at first than on the left side. This is due to the structure of the shutter, as can be seen in Appendix B. The overall opening speed on the other hand is not affected by this and seems to be linear with the offset.

2.4.1. Fast kinetics mode

The fast kinetics mode allows the image acquisition timings to be only dependent on the vertical shift speed, reducing the acquisition time significantly. In this mode, only a portion of the CCD is illuminated, while the dark parts of the chip will be used as a storage. This means, that as soon as an image is taken, the illuminated pixels are shifted vertically behind a mask, such that no photons can reach them anymore. When the chip is full or the user has finished their acquisition, the readout process is started.

As explained earlier in Chapter 2.2.1, the readout consists of first shifting a row into the readout register and then horizontally shifting them into the ADC. This is a very time con-

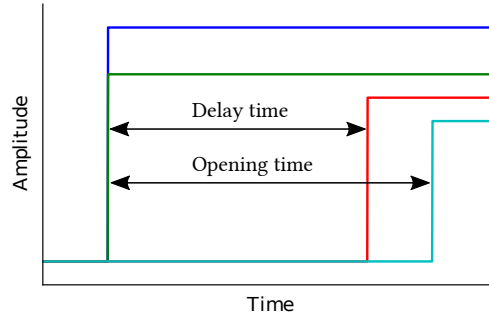


Figure 2.11.: **Shutter opening timings.** In the figure, the four signals show the timings until the shutter is open. The blue signal is the trigger that is sent from the user to start the opening sequence. The electronics add some minor delay, which is about a few nanoseconds until they send out the green trigger. After a certain delay time, which is most likely due to the coil needing to accumulate charge first, the red signal indicates when the shutter starts to open at the centre. The opening sequence is finished, when the shutter is fully open, which is shown by the turquoise line.

suming process, since the total readout time is described by

$$t_{ij} = iv_{ss} + (i - 1)j_{max}h_{ss} + jh_{ss}, \quad (2.6)$$

where v_{ss} and h_{ss} indicate the vertical and horizontal shift speed respectively. With this equation, the readout time until a pixel with the coordinates i and j is shifted into the readout register can be calculated. The readout of all pixels beforehand is also taken into account by j_{max} which is the width of the chip in pixels (1024 for the iKon M camera). As can be seen, in terms of the horizontal speed, the position of the last pixel in a row has a quadratic dependency. This is the dominant contribution to the readout time which is necessary to shift the pixels into the ADC. At this point, the speed should not matter too much anymore, because the experiment is already finished and the pixels are not illuminated anymore.

To set up the fast kinetics mode, there are several parameters that need to be set in advance.

- **Series length.** The number of images acquired before the readout phase is initiated.
- **Exposed rows.** The height of an image in the fast kinetic series.
- **Offset from bottom.** Number of rows from the bottom of the chip which are used as a temporary storage for the illuminated pixels.

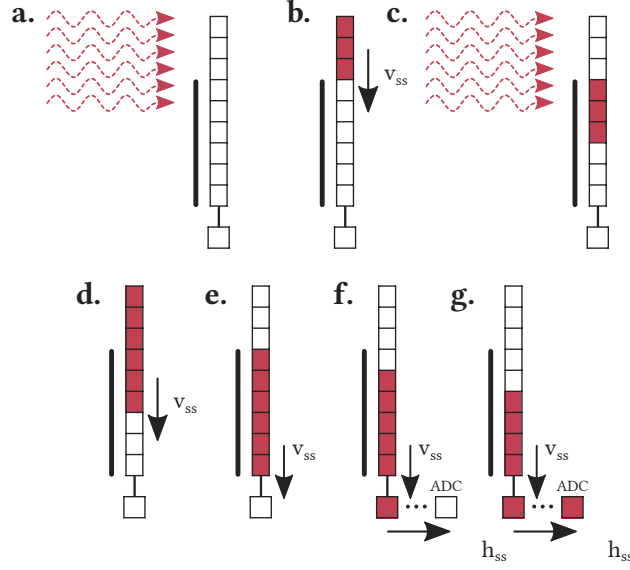


Figure 2.12.: **Schematics of imaging in the fast kinetics mode.** **a.** Laser photons excite electrons in pixels, creating the first absorption image. **b.** The illuminated pixels are shifted down behind a cover, while the laser is shut off. **c.** The second absorption image is taken without affecting the previous measurement. **d.** Both images are shifted down before starting the readout process, to not falsify the data with stray light. **e.** The readout process starts, the first row is shifted into the readout register. **f./g.** The first pixels are shifted into the ADC. It is repeated until the complete chip has been read out.

It can be easily seen that the parameters are limited by the height of the CCD in pixels:

$$\text{Exposed rows} + \text{Offset from bottom} \leq 1024 \text{ px} \quad (2.7)$$

and

$$\text{Series length} * \text{Exposed rows} \leq \text{Offset from bottom}. \quad (2.8)$$

Therefore, in our application, we chose an exposure of 204 px with an offset of 820 px, which gives a series length of 4. The pixels that are used as temporary storage of accumulated charge until the readout process begins need to be shielded from light. This is achieved with a slit in front of the camera that blocks part of the light beam implying on the CCD detector.

Since we will be imaging two atomic species, an absorption image for each species will be followed by two division images, as well as background images to subtract the noise. The fast kinetic mode now allows us to do this quickly, which is important when imaging multiple species together.

2.4.2. Frequency response of a slit

The slit is an optical element in the path, which will introduce diffraction — an effect due to the wave nature of light. Figure 2.13 sketches the systematics behind a plane wave approaching a slit.

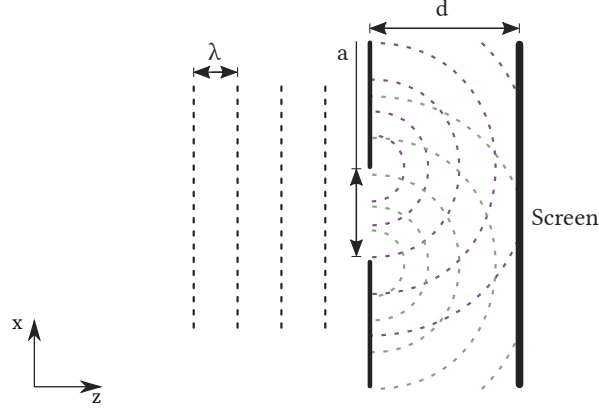


Figure 2.13.: **Diffraction on a single slit.** A planar wave with the wavelength λ approaches a slit with a width a , where the maxima of the wave are drawn as dashed lines. After the slit, Huygens principle is visualized by spherical waves, (purple and green), interfering with each other. The new wavefront is created where the waves intersect each other and is as such visible on the screen to the right.

The incoming planar wave can be described by the electric field:

$$E(z) = E_0 e^{-ikz}, \quad (2.9)$$

where the wave is propagating in the z -direction with the wave number k and an amplitude E_0 .

It is known from Huygens principle, that each point of a planar wave can be seen as the origin of a spherical wave. The spherical waves will interfere with each other and because of the superposition principle a new wavefront will be built up. The mathematical formalism of Huygens principle is simply the sum of all spherical waves, which for infinitely small distances is given as [14]:

$$E(x, z) = E_0 C \int_{\text{Slit}} \frac{1}{r^2} \exp(-ikr) dx'. \quad (2.10)$$

Here, $r = \sqrt{(x - x')^2 + z^2}$ is the radius of a two dimensional wave, C is a normalization constant. Taylor expansion of the radius for large z allows us, to substitute r^2 with z^2 and r

with $\frac{(x-x')^2}{2z} + z$, while also the approximation $(x - x') \ll z$ needs to be fulfilled. Therefore we get:

$$E(x, z) = E_0 e^{-ikz} \sqrt{\frac{ik}{2\pi z}} \int_{Slit} \exp\left(-\frac{ik}{2z}(x - x')^2\right) dx'. \quad (2.11)$$

The normalization was derived from the gaussian normal distribution $\exp(-\frac{ik}{2z}(x - x')^2)$. The integral will run over the slit size, with the origin in its middle and width a . We want to also write the exponential as a function of $-\frac{i\pi t^2}{2}$, due to the definition of fresnel integrals. The substitution follows:

$$\frac{k}{2z}(x - x')^2 = \frac{\pi t^2}{2}, \quad (2.12)$$

$$\Rightarrow t(x') = \sqrt{\frac{k}{z\pi}}(x - x'), \quad (2.13)$$

$$\Rightarrow dx' = -dt \sqrt{\frac{z\pi}{k}}, \quad (2.14)$$

such that the field becomes

$$E(x, z) = -E_0 e^{-ikz} \sqrt{\frac{i}{2}} \int_{t(-a/2)}^{t(a/2)} \exp\left(-\frac{i\pi t^2}{2}\right) dt. \quad (2.15)$$

The fresnel integrals are defined by

$$C(x) = \int_0^x \cos\left(\frac{\pi t^2}{2}\right) dt, \quad (2.16)$$

$$S(x) = \int_0^x \sin\left(\frac{\pi t^2}{2}\right) dt, \quad (2.17)$$

so that in combination with Eulers equations and splitting up the integrals the equation for the electric field is:

$$E(x, z) = E_0 \sqrt{\frac{i}{2}} e^{-ikz} [C(-a/2) - C(a/2) - iS(-a/2) + iS(a/2)]. \quad (2.18)$$

In a real experiment, the light on a CCD chip is detected as intensity, which is given by

$$I(x, z) = \frac{2|E(x, z)|^2}{\epsilon_0 c}, \quad (2.19)$$

with the dielectric constant ϵ_0 and the speed of light c .

To verify the theory, we set up a simple experiment with a collimated laser beam pointing at the CCD of the Andor camera, with a slit in between. The predicted interference fringes were observed in our experiment. Figure 2.14 shows the diffraction pattern that was recorded with the Andor camera and a slit. A large, collimated beam was used to illuminate a slit in front of the CCD detector, which casted the diffraction pattern on the chip. Using the theoretical description, we were able to extract the slit width $a = (2.470 \pm 0.013)$ mm and the distance to the chip $d = (11.0 \pm 0.3)$ mm.

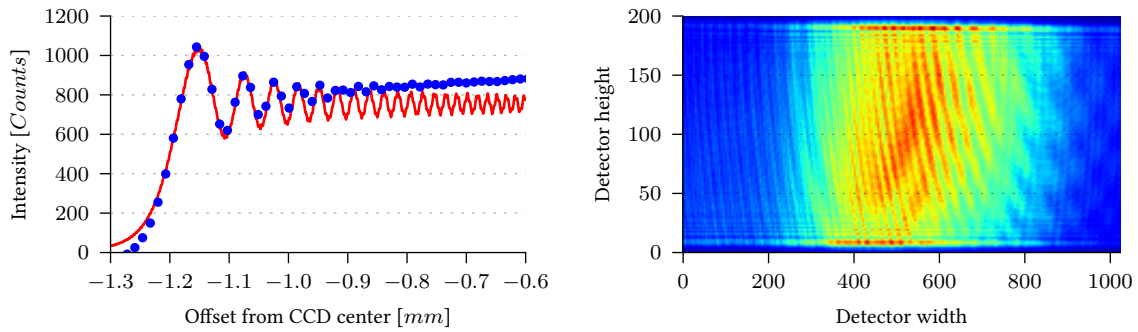


Figure 2.14.: **Measuring diffraction on a slit.** In order to characterize the diffraction by the slit, the CCD detector was placed as close as possible. The parameters were measured using a ruler and yielded distance $d = (10.9 \pm 0.5)$ mm, opening $a = (2.5 \pm 0.5)$ mm. The wavelength was $\lambda = 852$ nm found from the laser specifications. The blue curve is the experimental data, while the red curve was fitted, leaving the distance and opening as free fitting parameters. They were found to be $d' = (11.0 \pm 0.3)$ mm and $a' = (2.470 \pm 0.013)$ mm, which is in close agreement with the measured values. The residual deviation in the amplitude are caused by Gaussian intensity distribution of the laser beam used for illumination.

The experimental diffraction pattern matches the theory very well, since the deviations are in the expected regimes. The function washes out as it approaches the centre of the chip. This is due to the nature of the pixels, which only have a finite size and the fact, that the frequencies of the oscillations are lower on the outer ends. All oscillations that fit into one pixel are averaged, therefore diffraction is not visible.

This result helps to minimize the diffraction on the chip and serves as a starting point in order to optimize this effect.

2.4.3. Optimization of the masking setup

An important issue in the optimal placement of the slit is the distance from the CCD detector at which it is going to be placed. The frequency of the diffraction pattern depends on the distance d ; it is larger for smaller d , and smaller for larger d . In the limit $\frac{a^2}{d\lambda} \ll 1$ the diffraction pattern of a point source is recovered (Fraunhofer diffraction in the far field).

This dependence on the distance was also tested experimentally and the results are shown in Figure 2.15. The signal in one pixel is the average over all oscillations, that fit into the width of the pixel. Therefore the optimal slit position is close to the chip, since then the frequencies are maximal.

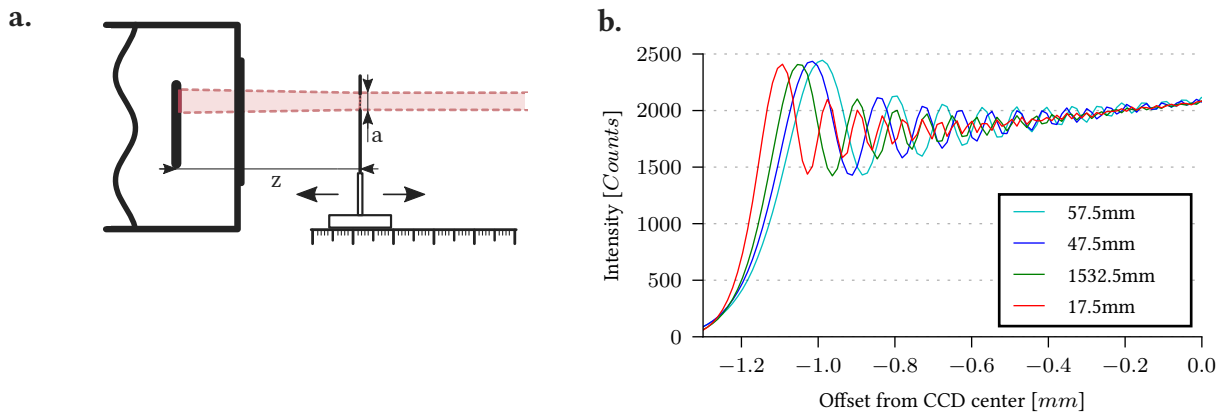


Figure 2.15.: **Distance dependent diffraction.** **a.** A slit was placed on a movable platform and diffraction was measured for various offsets z , while the slit opening a was kept constant. **b.** The diffraction frequency rises as the distance gets closer to the CCD.

3. Thermometry of an ultracold ideal fermi gas

The purpose of the camera is to measure scientifically important data from dense atomic clouds consisting of Lithium and Caesium. The improvement of the resolution in the whole imaging setup now allows to explore new attributes that could not be measured before and as an example, ultracold ideal Fermi Gases were chosen. They are apparent at very low temperatures, where only few atoms due to evaporative cooling are present. Their density distribution differs slightly from a gaussian form and from that, one can extract temperatures and atom numbers.

In order to image the atoms, a technique called absorption imaging is used which is explained in the next section, before the introduction to ideal Fermi Gases.

3.1. Absorption imaging

In order to find microscopic attributes of atoms, or systems of atoms, it is necessary to look at the atoms themselves. This is commonly accomplished using either fluorescence or absorption imaging [13]. In both cases, a laser beam is pointed at an atomic cloud, that is cooled and confined in a trap. In fluorescence imaging, the scattered light is collected, typically in a direction that is different than the illuminating beam. The intensity from the light through this method is not very high, since it is radiated in all directions. Therefore, long exposure times are required during which atoms can move and the information about the initial density and energy distribution is lost. Nevertheless, it is useful for single- and few-atom detection.

In contrast, in absorption imaging [15], the transmitted intensity of the imaging beam is recorded. Without atoms, one would see a beam profile of the laser beam. With atoms, a shadow is visible due to the atoms "blocking" the light. This is accomplished, by correctly tuning the laser to a resonance frequency of the atoms, which enables them to absorb the

cite phd thesis of Serwane, Selim, Science 332, 336 (2011)

light, exciting them to a higher state. Through spontaneous emission, the atoms will decay, making it possible to excite them once again. This method works well, when the "signal" from the absorbed light is significantly larger compared to the noise sources.

There are a set of optical elements in the imaging path, like lenses to collimate the image and refocus it, or mirrors to guide the light into the camera. Since the surfaces will most likely introduce errors into the imaging, for example from impurities or dust, only the absorption image will not suffice to gain reliable data. This is compensated by taking a total of three pictures, in order to extract only the relevant information from the image.

This can be understood when looking at the light intensity I_{CCD} reaching the camera. The atom cloud has an optical density OD , therefore the intensity can be written as [13]

$$I_{CCD} = I_0 e^{-OD} + I_{back}, \quad (3.1)$$

where it decreases from the incident laser intensity I_0 due to light scattering by atoms. The intensity I_{back} describes the background signal, that is found when the CCD is not being illuminated by a laser such as readout noise, dark noise or stray photon light. All the interesting attributes of atoms are found by looking at the optical density, therefore in order to extract that, a background frame is subtracted from the absorption image and the laser profile divided, leaving

$$\frac{I_{CCD} - I_{back}}{I_0} = e^{-OD}. \quad (3.2)$$

The laser intensity I_0 is measured in a separate frame, containing the laser intensity $I'_0 = I_0 + I_{back}$ and also the background I_{back} . Finally, the equation yields

$$\frac{I_{CCD} - I_{back}}{I'_0 - I_{back}} = e^{-OD}. \quad (3.3)$$

three images
as example

From the resulting optical density, one can now conclude, for example, atom density distributions, atom numbers or excitation rates.

3.2. Density distributions of ideal Fermi gases

Ideal Fermi gases offer a new aggregate, that is complementary to Bose-Einstein condensate. They are found on the Bardeen-Cooper-Schrieffer side of the Feshbach resonances as a

polarized species with only one spin component. At very cold temperatures, they start to differ from a gaussian distribution which is further investigated in this chapter.

The distribution of the atoms depends on the fraction of their temperature to the Fermi temperature [16] T_F . For $T/T_F \gg 1$, the atoms will follow a gaussian distribution, which can be identified using the gaussian radius:

$$\sigma_i = \sqrt{\frac{2k_B T}{m w_i^2}}, \quad (3.4)$$

with the mass of Lithium m and the trapping frequency w_i in the direction of the radius. Due to the alignment of the dipole trap, the cloud will not have a spherical shape and will therefore have different radii.

In the degenerate regime however, for $T/T_F \ll 1$, the radius is described by the Fermi radius

$$R_{Fi} = \sqrt{\frac{2E_F}{m w_i^2}}, \quad (3.5)$$

using the Fermi energy E_F , due to the fermions filling up the eigenstates of the potential.

It is therefore suggested [16] to use a unified radius, as the temperatures are not known a priori:

$$R_i^2 = \frac{2k_B T}{m w_i^2} f(e^{\frac{\mu}{k_B T}}). \quad (3.6)$$

The interpolation function $f(x)$ is hereby:

$$f(x) = \frac{Li_1(-x)}{Li_0(-x)} \quad (3.7)$$

where Li_n is polylogarithm and can be defined as

$$Li_s(z) = \sum_{k=1}^{\infty} \frac{z^k}{k^s}. \quad (3.8)$$

In our case, we integrate over all but one axes, and therefore find the fitting function for the atom numbers:

$$n_{1D}(x) = n_{1D,0} \frac{Li_{5/2} \left(\pm \exp \left[q - \frac{x^2}{R_x^2} f(e^q) \right] \right)}{Li_{5/2}(\pm e^q)}. \quad (3.9)$$

The derivation can be found in [16]. The parameter $q = \frac{\mu}{k_B T}$ can be extracted, which contains information about the chemical potential μ and the temperature T .

This parameter can then be used to calculate the degeneracy parameter:

$$\frac{T}{T_F} = [-6Li_3(-e^q)]^{-1/3}. \quad (3.10)$$

To compensate for the finite resolution of the chip on the camera, the cloud can be imaged after a short time of flight t . The temperature can then be calculated from the dynamics as

$$k_B T = \frac{1}{2} m w_i^2 \frac{R_i^2}{1 + w_i^2 t^2} \frac{1}{f(e^q)}. \quad (3.11)$$

3.3. Finding properties of the Fermi gas

As seen before, from an ideal Fermi gas, one can deduce the temperature and Fermi temperature of a gas. In order to implement ideal Fermi gases, Lithium was prepared in an optical dipole trap. The power was ramped down several times at the Feshbach resonance $B = 896$ G, until only a fraction of the atoms remained. The spin-down atoms received a short laser pulse, so that the cloud only consisted of spin-up ^6Li .

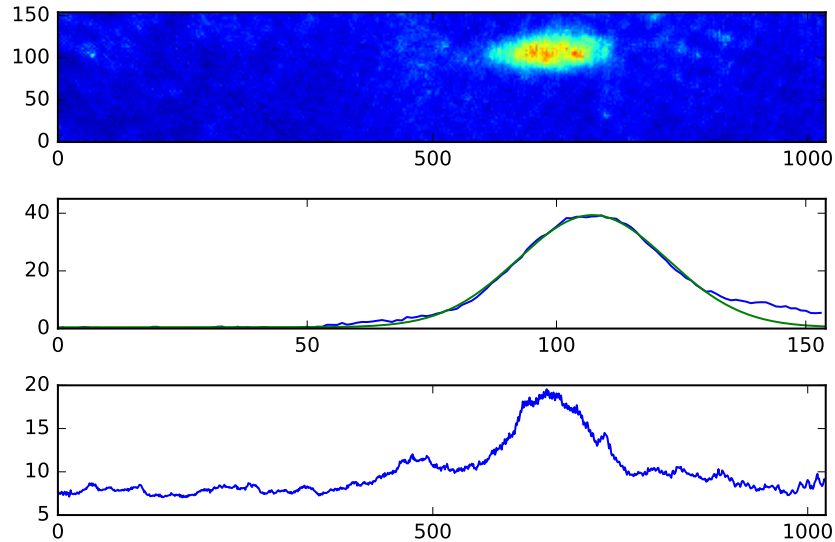


Figure 3.1.: **saywhat.** plt pls.

The measurement was executed for three different powers of the laser beam. The atoms were imaged in situ (in the trap) and after a short time of flight of 1 ms. In order to fit the

function Equation 3.9, the absorption image yielding the optical density as seen in Figure 3.1 was integrated over one axis, giving a one dimensional distribution for the atoms.

As the imaging is not properly optimized at this point, the distribution in the horizontal (x-axis) direction does not represent the theory well, even after averaging over several acquisitions. This problem could not be resolved due to timing constrictions, therefore the fit parameters were extracted solely from the vertical direction.

On the vertical axis, a common systematic error was a small peak for high values, which did not correspond to the atoms and is probably due to the lenses, which are not properly aligned. This was excluded from the fit. The results are found in [reftable].

reftable

4. Conclusion and outlook

Using the new setup, it is now possible to take images with low noise even at few atom numbers. With the high quantum efficiency, atoms are detected in more detail than before. The water cooling helps to cool the chip of the camera uniformly to allow minimal dark noise, which then separates the signal from the atoms from the background even at low atom numbers. Using the fast kinetics mode, the next step would be to take images of Lithium and Caesium consecutively. Previously, images of both species had to be taken from different clouds or axes.

The polarized Fermi gas already shows the potential of the new system. Building on that, one could now as well create a Fermionic superfluid from Lithium, using a spin imbalance. This new aggregate has a unique structure, where there is an unpolarized core surrounded by the majority spin component. The latter is then as well described by the same formulas as the polarized Fermi gas, which was also discussed in this thesis.

cite some
superfluid
paper

A. Acquisition sequence

During absorption imaging, the Andor camera will take eight images, which consist of two absorption, two division and four background images (or if the user specifies to only image one species, then there is one absorption and division image and six background images).

As the camera is only interested in the rising edge of the signal, the trigger length is not really important. But the chip will be exposed for the duration of the exposure time, after which the laser should be turned off. The, in our case, 204 pixels are then shifted downwards, taking t_{vshift} , which can be calculated by multiplying the vertical shift speed with the pixel height. For the fastest shift speed, this would result in

$$t_{vshift} = 2 \mu\text{spx}^{-1} * 204 \text{ px} = 408 \mu\text{s}, \quad (\text{A.1})$$

for the slowest speed we find

$$t_{vshift} = 64 \mu\text{spx}^{-1} * 204 \text{ px} = 13.1 \text{ ms} \quad (\text{A.2})$$

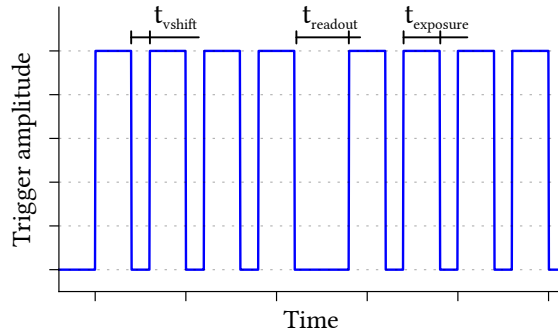


Figure A.1.: **Trigger signals.** The triggers are sent as a rectangular signal. The camera will interpret the rising edge of the signal as the starting point of the exposure. Here, t_{vshift} is the time it takes to shift the illuminated pixels downwards, $t_{readout}$ is the time it takes the chip to read out the data and $t_{exposure}$ the exposure time.

Trigger number	Trigger time
1	0
2	t_{trig}
3	$2 * t_{trig}$
4	$3 * t_{trig}$
5	$4 * t_{trig} + t_{readout}$
6	$5 * t_{trig} + t_{readout}$
7	$6 * t_{trig} + t_{readout}$
8	$7 * t_{trig} + t_{readout}$

Table A.1.: **Trigger timing.** The trigger signals are limited by the readout speeds of the chip. This table lists the minimal timings necessary to fully read out the chip. Between signal 4 and 5, the illuminated pixels are still shifted downwards, before the chip is read out. This prevents introducing unnecessary errors, as all relevant pixels are then behind the cover of the slit.

Therefore the four subsequent signals can be taken after

$$t_{trig} = t_{vshift} + t_{exposure}. \quad (\text{A.3})$$

Four images are acquired until the readout process is started. This takes significantly longer, but as explained in Chapter 2.4.1 can be calculated via:

$$t_{ij} = i * v_{speed} + (i - 1) * j_{max} * h_{speed} + j * h_{speed}. \quad (\text{A.4})$$

For the whole chip, this gives us $i = 4 * 204 = 816$, $j = 1024$, $j_{max} = 1024$, therefore

$$t_{readout} = 816 * v_{speed} + 816 * 1024 * h_{speed}. \quad (\text{A.5})$$

The trigger signals can finally be calculated via Table A.1.

B. Design of the shutter

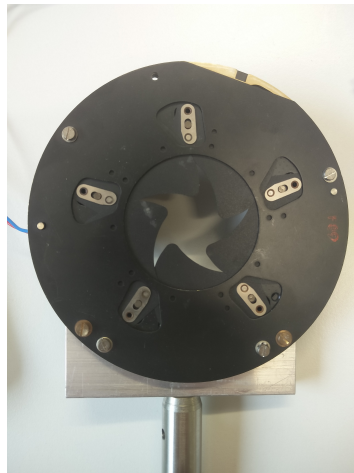


Figure B.1.: **Front view of the shutter.** Five fans are visible on the shutter as well as the mechanical guides, used to pull the fans outwards.



Figure B.2.: **Back view of the shutter.** The magnet is connected to the lever, which pulls the fans outwards.

C. Shutter circuit

The circuit given on the next page is the complete circuit as opposed to Figure 2.7. The key difference here, is that the actual board contains two input/output ports and there are also converters (visible in the top-left corner), that keep the voltage at 24 V, which is necessary in order to use the MOSFET driver.

The breadboard on the bottom are built into the circuit, so that custom additions can be made, although this is not used in our case. This circuit is commonly used in various labs in the institute, if high currents are required and was designed by Gerhard Zuern.

Reference
Gerhard

D. Setup of the custom slit

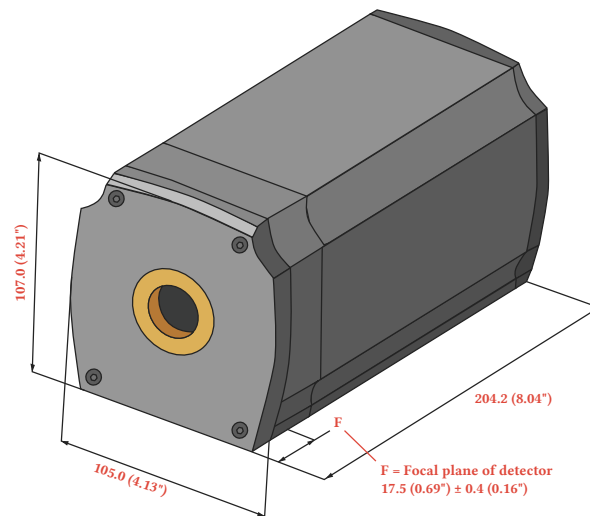


Figure D.1.: **Drawing of the camera.** As can be seen in this drawing, the CCD chip is first hidden behind a cover, that also includes an internal shutter and then offset by an additional 5 mm.

The cover, which can be seen in Figure D.1 has a width of 12.5 mm, which adds additional space before the chip. The cover is mainly for a manual cap to cover the chip, when the camera is not used, and an internal shutter. Since we knew, that the internal shutter was not needed, we were able to remove the cover bringing the slit closer to the chip. Images of under the cover can be found in Appendix .

A/B/C? Im-
age?

The holes for M4 screws were already there, so a custom plate was built on which the slit could be mounted on. In the technical drawing in Figure ??, the centre-most holes are reserved for the slit, which can be moved up and down to select the appropriate height needed for the imaging. The plate also gives the opportunity to move the whole set with the long holes in the outer-most edges.

It is also important to note, that the camera is very sensitive to stray light due to the high quantum efficiency. Therefore it is necessary to cover the laser path. The plate offers an-

other set of screw holes, which will hold a SM2-mount, therefore eliminating any gap that could allow photons to reach the camera externally.

The long path of SM2 tubes then only allows the smallest amount of stray light to enter the camera, which will have significantly lower intensity than the actual absorption image.

List of Figures

2.1. Imaging path	3
2.2. Schematic design of a CCD array and its readout	5
2.3. Shifting charges in a CCD detector	5
2.4. Dark noise accumulation on the chip	8
2.5. Gradient on the chip during readout	8
2.6. Readout noise	9
2.7. Electronic circuit to control the shutter	12
2.8. Probing the shutter for dynamics	13
2.9. Shutter characterization	14
2.10. Sample dynamics	15
2.11. Shutter opening timings	16
2.12. Schematics of imaging in the fast kinetics mode	17
2.13. Diffraction on a single slit	18
2.14. Measuring diffraction on a slit	20
2.15. Distance dependent diffraction	21
3.1. saywhat	26
A.1. Trigger signals	31
B.1. Front view of the shutter	33
B.2. Back view of the shutter	33
D.1. Drawing of the camera	37

Bibliography

- [1] University of Heidelberg. <http://dc.zah.uni-heidelberg.de/lswscans/res/positions/q/form>. Accessed: 28.03.2016.
- [2] R. Pires, J. Ulmanis, S. Häfner, M. Repp, A. Arias, E. D. Kuhnle, and M. Weidemüller. Observation of Efimov resonances in a mixture with extreme mass imbalance. *Phys. Rev. Lett.*, 112:250404, Jun 2014.
- [3] J Ulmanis, S Häfner, R Pires, E D Kuhnle, M Weidemüller, and E Tiemann. Universality of weakly bound dimers and efimov trimers close to Li–Cs Feshbach resonances. *New Journal of Physics*, 17(5):055009, 2015.
- [4] J. Ulmanis, S. Häfner, R. Pires, F. Werner, D. S. Petrov, E. D. Kuhnle, and M. Weidemüller. Universal three-body recombination and Efimov resonances in an ultracold Li-Cs mixture. *Phys. Rev. A*, 93:022707, Feb 2016.
- [5] Marc Repp. *Interspecies Feshbach Resonances in an Ultracold, Optically Trapped Bose-Fermi Mixture of Cesium and Lithium*. PhD thesis, Ruperto-Carola-University of Heidelberg, 2013.
- [6] Andor. *Andor iKon M manual*.
- [7] Chris McFee. An introduction to ccd operation. http://www.mssl.ucl.ac.uk/www_detector/ccdgroup/opttheory/ccdoperation.html. Accessed: 23 March 2016.
- [8] Allied Vision. Guppy f-038 datasheet.
- [9] M. Repp, R. Pires, J. Ulmanis, R. Heck, E. D. Kuhnle, M. Weidemüller, and E. Tiemann. Observation of interspecies ${}^6\text{Li}$ - ${}^{133}\text{Cs}$ Feshbach resonances. *Phys. Rev. A*, 87:010701, Jan 2013.
- [10] Carmen Renner. Design and characterization of a dual-wavelength high-resolution imaging system. Master’s thesis, University of Heidelberg, 2014.

- [11] *Skript zum Versuch F36 – Teil I "Wellenfrontanalyse mit einem Shack-Hartmann-Sensor" des Fortgeschrittenen-Praktikums II der Universität Heidelberg für Physiker.*
- [12] Timo Bastian Ottenstein. A new objective for high resolution imaging of Bose-Einstein condensates. Master's thesis, University of Heidelberg, 2006.
- [13] Simon Murmann. Toward the exact detection of mesoscopic atom numbers. Master's thesis, University of Heidelberg, 2011.
- [14] Daniel A. Steck. Quantum and atom optics. Available online at <http://steck.us/teaching>, April 2014. revision 0.9.10.
- [15] Stephan Helmrich. Improving optical resolution by noise correlation analysis. Master's thesis, 2013.
- [16] Wolfgang Ketterle and Martin W. Zwierlein. Making, probing and understanding ultracold fermi gases. In W. Ketterle M. Inguscio and C. Salomon, editors, *Proceedings of the International School of Physics Enrico Fermi*, volume Course CLXIV, page 95. IOS Press, Amsterdam, 2008.